



medi-inn PRO

Urine bags


Closed urinary drainage system



Product Benefits

- ✓ Closed urinary drainage system with needle-free urine sampling port
- ✓ Universal stepped connector
 - ensures secure catheter connection
- ✓ Needle-free urine sampling port helps prevent needlestick injuries and allows easy cleaning thanks to its flush design
- ✓ Large drip chamber with integrated anti-reflux valve and constant pressure equalisation through the hydrophobic-coated microfilter reduces the risk of urinary tract infections
- ✓ Sturdy housing ("doghouse") ensures secure fixation of the drainage outlet
- ✓ Hygienic, convenient one-handed operation thanks to urine drainage with cross valve
- ✓ Secure fixation with
 - flexible, stable universal double hanger
 - additional strap with hanging loop
- ✓ Kink-resistant, flexible tubing with clamp and sheet clip

Product Details

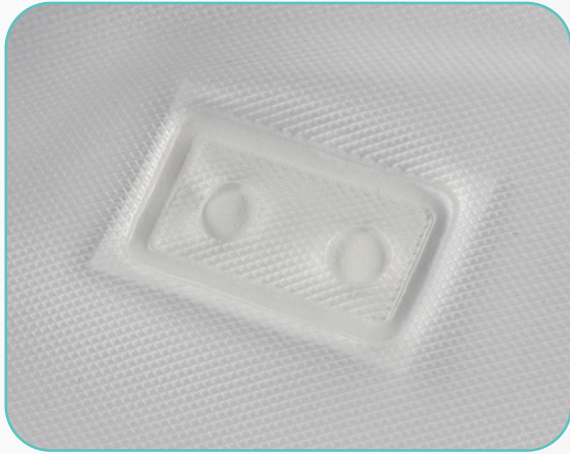

Liquid-tight, hydrophobic-coated microfilter



Drip chamber with integrated anti-reflux valve and hydrophobic-coated microfilter



Connector with needle-free urine sampling port helps prevent needle-stick injuries



Urine drainage with cross valve for hygienic and convenient one-handed operation

Ordering Information

Item No.	Product	Tube Length	Type	Drain Valve	Contents per Shipping Unit (PU)
N19008	medi-inn PRO Urine Bag 2000 ml	120 cm	Needle-free	Cross valve, vertical	25 pieces

✓ Class Is sterile medical device in accordance with Regulation (EU) 2017/745

✓ Latex-free

✓ Shelf life: 5 years



WWW.BODY-PRODUCTS.DE

Your Contact:

BODY PRODUCTS
Relax Pharma und Kosmetik GmbH

Alfred-Nobel-Straße 1-3
50226 Frechen

Phone:	+49 2234 91165-0
Fax:	+49 2234 91165-23
Mail:	info@body-products.de