



medi-inn PRO

Urine Measurement System

Closed single-chamber urine measurement system





Product Benefits

- ✓ **Closed single-chamber urine measuring system with fine graduation on the urimeter for precise monitoring and measurement of the urine volume**
 - Finely graduated urimeter for precise urine output measurement
 - Additional drainage outlet on the urimeter
- ✓ **Universal stepped connector**
 - ensures secure catheter connection
- ✓ **Large drip chamber with integrated anti-reflux valve and constant pressure equalisation through the hydrophobic-coated microfilter reduces the risk of urinary tract infections**
- ✓ **Sturdy housing ("doghouse") ensures secure fixation of the drainage outlet**
- ✓ **Needle-free sampling port**
 - reduces the risk of needlestick injuries
 - allows easy cleaning due to its smooth design
- ✓ **Large measurement scale with clearly readable markings**
- ✓ **Hygienic, convenient one-handed operation thanks to urine drainage with cross valve**
- ✓ **Secure fixation with**
 - flexible, stable universal double hanger
 - additional strap with hanging loop
- ✓ **Kink-resistant, flexible tubing with locking clamp and sheet clip**
 - Anti-kink spiral at the proximal tube end reduces the risk of kinking


Product Details



Liquid-tight, hydrophobic-coated microfilter



Drip chamber with integrated anti-reflux valve and hydrophobic-coated microfilter

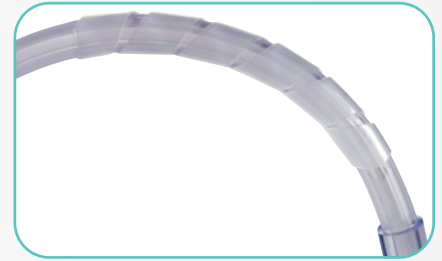

Anti-kink spiral at the proximal tube end reduces the risk of kinking



Connector with needle-free urine sampling port helps prevent needlestick injuries



Urine drainage with cross valve for hygienic and convenient one-handed operation

Ordering Information

Item No.	Product	Collection Bag Capacity	Urimeter Capacity	Tube Length	Type	Drain Valve	Contents per Shipping Unit (PU)
N18984	medi-inn PRO Urine Measurement System	2.600 ml	500 ml	150 cm	Needle-free	Cross valve, vertical	10 pieces

✓ Class Is & Im sterile medical device with measuring function in accordance with Regulation (EU) 2021/745

✓ Latex-free

✓ Shelf life: 5 years



WWW.BODY-PRODUCTS.DE

Your contact:

BODY PRODUCTS
Relax Pharma und Kosmetik GmbH

Alfred-Nobel-Straße 1-3
50226 Frechen

Phone: +49 2234 91165-0
Fax: +49 2234 91165-23
Email: info@body-products.de